SOEHNLE



BODY CONTROL CONTOUR F3



BEDIENUNGSANLEITUNG
OPERATING INSTRUCTIONS
MODE D'EMPLOI
ISTRUZIONI PER L'USO
GEBRUIKSAANWIJZING
INSTRUCCIONES DE MANEJO
MANUAL DE INSTRUÇÕES
BRUKSANVISNING
BRUGSANVISNING
KÄYTTÖOHJEET
KEZELESI UMUTATO
INSTRUKCJA OBSŁUGI
NÁVOD K POUŽITÍ
ИНСТРУКЦИИ ПО ЭКСПЛУАТАЦИИ
KULLAMA KILAVUZU
Οδηγίες λειτουργίας
NAVODILA ZA UPORABO
UPUTA ZA UPORABU



1.

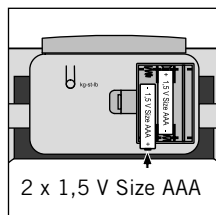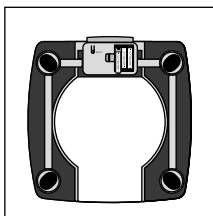2.

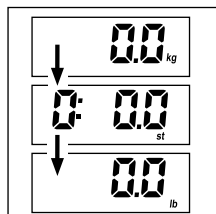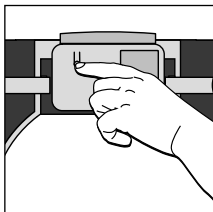3.

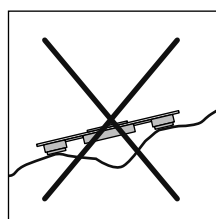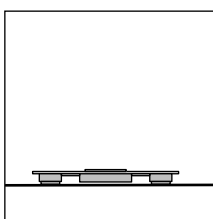D

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General information

Guarantee

Soehnle guarantees to rectify defects caused by material or production faults for 3 years from date of purchase free of charge, by means of repair or exchange. Please ask retailer to fill in and stamp guarantee slip (on page 23) at time of purchase. When making a claim under guarantee, please return scale to your retailer with guarantee slip.

Conformity

This appliance is radio interference-free in accordance with the applicable EC directive 2004/108/EC. You will find the EU compliance declaration on www.soehnle.com.

Although the appliance has been tested and its emissions do not exceed the permitted radio frequency emission levels as laid down in the above directives, there is no guarantee that interference will not occur in certain situations. Possible interference can be detected by switching the appliance on and off at the same time as checking the radio or TV reception. You can try to remove interference by carrying out one or more of the following steps:

- Point the appliance in question or its antenna in a different direction.
- Increase the distance between the appliance in question and the scale.
- If necessary, call in an experienced radio/TV engineer.

IMPORTANT! Making modifications to this product's electronics or housing requires express consent from Soehnle; ignoring this regulation may result in the user's permission to operate the appliance being withdrawn.

Electromagnetic influences

Under extreme electromagnetic influences, e.g. operating a radio-controlled appliance in the immediate proximity of the appliance, this may cause a detrimental effect to the display value. When the cause of the disturbance ceases, the appliance can be used again as intended, or it may be necessary to switch it off and on or temporarily remove the batteries.

Statutory notice for disposing batteries

Information: all batteries originally included with the product are non-hazardous! Batteries should not be disposed of with household waste. You have a legal duty as a consumer to return used batteries. You can hand in your old batteries at public collection points in your community or at any place selling batteries of the type in question.

Note: You will find these symbols on batteries containing hazardous materials:

Pb = battery contains lead

Cd = battery contains cadmium

Hg = battery contains mercury

Disposal of used electrical and electronic appliances

This symbol on a product or its packaging means that this product should not be treated as normal household waste, but must be handed in at a recycling collection point for electrical and electronic appliances. You can find out more information from your local council, the municipal waste companies or the shop where you purchased the product.

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LEIFHEIT AG